


Thank you for choosing our product !

The Ora Medical team is excited to accompany you as you take your first steps with Levity !

by ORA Medical Inc.

Share your experience with us on Instagram by tagging @ora_medical.

Unpacking the Levity

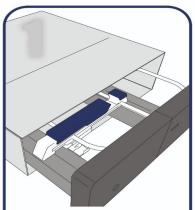


Take the device out of the box and make sure all items included in your order are present.

Tip: Lay the package on its widest side to slide the device out of the box.




Loosen the levers on the sides and unfold the upper frame.

Tip: Hold the bottom frame with one hand and pull the upper frame with your other hand.

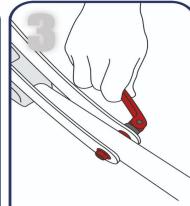

Ensure that the sliders are securely engaged in the safety grooves and tighten the levers on the sides.


Tip: Press on the sliders to ensure they are well engaged in the grooves.

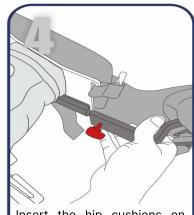

Insert the hip cushions on each side onto the harness rail and tighten the handles to secure their position.

Tip: Make sure the harness is untangled before installing it.

Installing the User

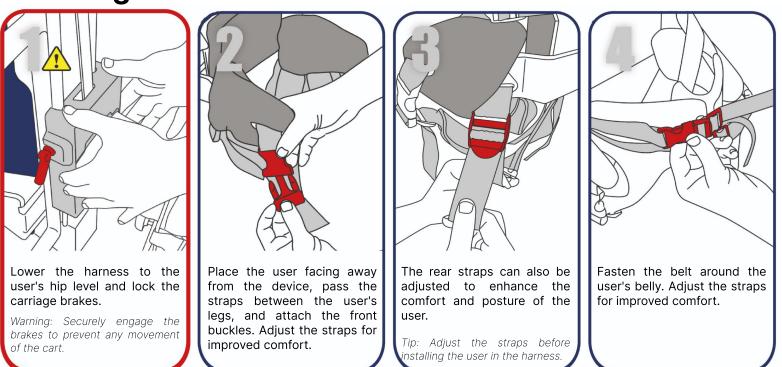

OrthoPediatrics Corp. is an exclusive distributor of Ora Medical, Inc.'s Levity Device "(Levity") worldwide (excluding Canada). OP Specialty Bracing products and products distributed by OrthoPediatrics Corp. and its subsidiaries should be used under the guidance of a qualified healthcare professional. Individual results may vary. Levity is intended to provide adjustable weight support for people with walking disabilities who lack strength, good balance, or endurance. Please consult your pediatrician or orthopedic specialist for professional advice, including product warnings, precautions, side effects and contraindications. This document is for informational purposes only and does not constitute medical advice. Always follow your doctor's recommendations and instructions.